
VIVIX-S Portable, Wired

Portable Flat Panel Detector
For Digital Radiography



VIVIX-S Portable, Wired is a Vieworks' flat panel digital radiography cassette system with 14" x 17" coverage area for general radiographic applications using its unique image processing system and proprietary flat panel detector. With the same size as a film or a CR cassette, it is an excellent solution for upgrading conventional X-ray system. User friendly imaging software, VXvue™ is compliant with DICOM 3.0 standard and provides consistent image quality at a dramatically reduced dose and faster image information with optimized algorithm for each different study. It only takes a few simple steps to acquire and transmit images to the DICOM server through Gigabit Ethernet.


Features

- * Wide active area of 14" x 17"
- * High spatial resolution with 140um pixel array
- * Stable and reliable automatic exposure detection (Anytime™)
- * Viewer software running on Windows™ OS (VXvue™)
- * Communication interface through Gigabit Ethernet (1000 BASE-T)
- * Simple and easy integration with all kinds of digital radiography system




Configuration
 Detector FXRD-1417SA/B
 System Control Unit FXRS-02A
 Designed and manufactured by Vieworks in Korea

Acquired Images

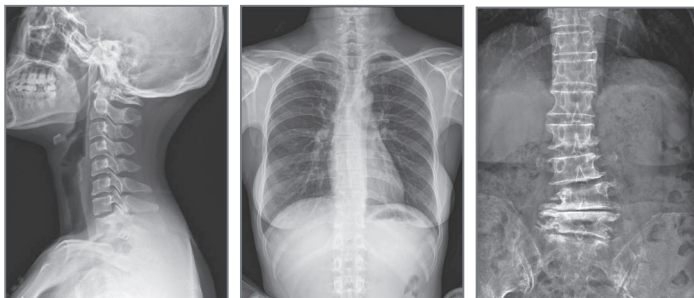

Drawing

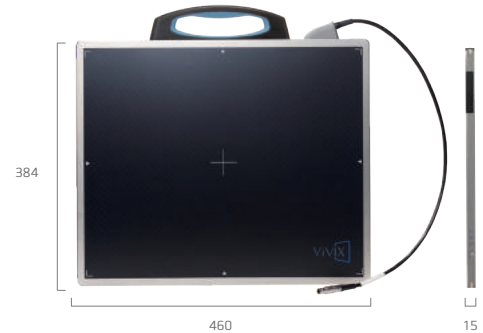

Technical Specification

Application	General radiography
Technology	Flat panel detector : a-Si TFT with PIN diode
Scintillator	CsI:Tl / Gd ₂ O ₂ S:Tb
Pixel Pitch	140um x 140um
Pixels	2,560 x 3,072 pixels
Image Size	14 x 17 inches (35 x 43cm)
A/D Conversion	14 bit
Grayscale	16,384 steps
X-ray Voltage Range	40 ~ 150kVp
X-ray Generator Interface	Line trigger : DR Trigger Mode Auto trigger : AED (Automatic Exposure Detection) Mode
Data Interface	Gigabit Ethernet (1000BASE-T) via PoE (Power over Ethernet)
Extension Cable	GigE communication and power supply (7m)
Dimensions	460 (W) x 384 (L) x 15 (T) mm
Weight	Approx. 2.9kg (GADOX) / 3.1kg (CsI)
Operating Environment	15 ~ 35°C, 30 ~ 85% RH (non-condensing)
Power	DC24V, 0.5A

* Specifications are subject to change without prior notification.

FDA Approved

1 Chapin Road, Unit 3, Pinebrook, New Jersey 07058
 ♦Tel: 800-456-7800 ♦International : 201-298-2980 ♦Fax: 888-437-9729
 Email: info@medlinkimaging.com

www.medlinkimaging.com